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EN DIGITAL THERMOMETER model Thermo Soft USER'S MANUAL

- Measuring time: 60±10 seconds - oral, rectal; 100±20 seconds - axillary
- Flexible tip
- Last measurement memory
- Fever alarm
- Auto shut-off
- Water proof

Thank you for choosing TM Gamma® Digital Thermometer Model Thermo Soft. We are sure that having appraised worthy the high quality and reliability of this device you will become a regular user of the products of English Trademark Gamma®. Before starting to use the thermometer please study the instructions carefully.

IMPORTANT SAFETY INSTRUCTIONS

1. The thermometer is suitable exclusively for measurement of the body temperature!
2. The minimum measurement time until the signaling tone «beep» must be maintained without exception!
3. The thermometer contains small parts (battery, etc.) that can be swallowed by children. Therefore never leave the thermometer unattended to children.
4. Do not store the unit where it will be exposed to direct sunlight, dust or humidity. Avoid extreme temperatures, dropping or subjecting your thermometer to strong shocks.
5. Do not attempt to disassemble the unit, except to replace the battery.
6. Avoid ambient temperatures above 55 °C. NEVER boil the thermometer!
7. Store the unit in the protective case when not in use.

 **Batteries and electronic instruments must be disposed of in accordance with locally applicable regulations, nor with domestic waste.**

FEATURES OF YOUR DIGITAL THERMOMETER

HOW TO TAKE YOUR TEMPERATURE

1. Press ON/OFF button to activate. The unit would beep and display «188.8 °C». This LCD display test will go on for around 2 seconds.
2. Then, the temperature taken last time would be displayed with a «M» mark in the right corner.

3. Then, 37.0°C ("Lo" and "Co" flashing) will appear.
4. Put the thermometer to places you want to take temperature, into mouth, under armpit or inside the anus.
5. The beep will sound when temperature measuring is complete (when taking an oral and rectal temperature). If the signalling tone sounds and «C» is no longer flashing, this means that the measured increase in temperature is now less than 0.1 C in 16 sec. When taking an axillary temperature, ignore the beep sound and wait 3-5 minutes else before taking the reading.
6. Press ON/OFF button to turn off the unit. (The unit will automatically power off after approximately nine minutes after use.)

MEASURING METHODS

| Measuring method | Position | Measuring time |
|---|---|-----------------|
| Taking an oral temperature (in the mouth). Position the thermometer in one of the two pockets under the tongue, to the left or right of the root of the tongue. The measuring sensor must be in good contact with the tissue. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/exhaled air. Do not eat or drink cold or hot products before measurement. |  | 60 ±10 seconds |
| Taking a rectal temperature (in the anus). This is the most reliable measuring method, and is especially appropriate for infants and small children. Apply a water-soluble lubricant to the probe cover. Gently insert the probe (MAXIMUM 1/2") into the rectum. |  | 60 ±10 seconds |
| Taking an axillary temperature (under the arm). From a medical viewpoint, this method will always provide inaccurate readings, and should not be used if precise measurements are required. |  | 100 ±20 seconds |

Note: Rectal temperatures are generally 0.5°C /1.0°F higher than oral readings. Axillary temperatures are generally 0.5°C /1.0°F lower than oral readings.

CLEANING INSTRUCTIONS

Before and after each use, clean the thermometer with a soft cloth and ethyl alcohol diluted with water, or wash with COLD soapy water. Do not soak for an extended period of time or sterilize by boiling method, gas or steam autoclaves.

REPLACING THE BATTERY

The battery in your digital thermometer needs to be replaced when «□» appears in the lower right corner of the display. Replace the battery with a 1.55V, LR41 or SR41. To replace the battery:

1. Pull the battery cover away from the unit.
 2. Use a non-metal pointed tool to remove the battery and insert the new battery with the positive(+)-side facing up.
 3. Carefully slide the battery cover back into place.
- Note:** Please properly replace the battery away from small children and heat.

SPECIFICATIONS


| | |
|------------------------------------|--|
| Measurement range | 32.0°C ~42.9°C |
| Low temperature display | Temperature <32.0°C: display: Lo°C |
| High temperature display | Temperature ≥42.9°C: display: Hi°C |
| Display resolution | 0.1°C |
| Accuracy | 34.0°C ~42.9°C ±0.1°C other range: ±0.2°C |
| Power Consumption | 0.15 milliwatts in measurement mode |
| Normal temperature alarm (≤37.8°C) | Ten quick beeps for approximately 10 seconds |
| Fever alarm (>37.8°C) | Three quick beeps followed by a pause for approximately 10 seconds |
| Auto power off | Approx. 9 minutes |
| Battery | LR41 or SR41 |
| Life of battery | More than 200 hours of continuous operation |
| Operating condition | Temperature: 10°C~40°C Humidity: 15%~95%RH |
| Storage condition | Temperature: -25°C~55°C Humidity: 15%~90%RH |

* The specifications and the outer design are subject to change for improvement reasons without previous notice.

WARRANTY

The digital thermometer is guaranteed for 3 years from date of purchase. The guarantee does not apply to damage caused by improper handling, accidents, not following the operating instructions or alterations made to the instrument by third parties. The guarantee is only valid upon presentation of the guarantee card filled out by the representative of the sales organization and its stamp and safety of the packaging.

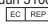Hangzhou Universal Electronic Co., Ltd. 38 Yangjiatang, Sandun 310030 Westlake District, Hangzhou China.  Prolinx GmbH, Brehmstr. 56, 40239 Duesseldorf, GERMANY.

 Prolinx GmbH
Brehmstr.56,40239 Duesseldorf, Germany

Gamma is a registered Trademark